

## Congratulations on your purchase of PediWow™


The most innovative safe and effective calluses, dry skin and corns remover EVER!

PediWow™ is a Professional grade Do-It-Yourself tool designed to quickly and easily remove calluses, dead skin and corns in any area and at the same time being safe and gentle.

It works well on your feet, toes, hands and fingers too!

PediWow™ (Fig. 1) even works great on both dry and wet skin. This device has two different exfoliating tips. The Power Tip (Fig. 1a) works on wide areas while the Precision Tip (Fig. 1c) works on small and difficult to reach areas.

Fig.1 PediWow™ housing

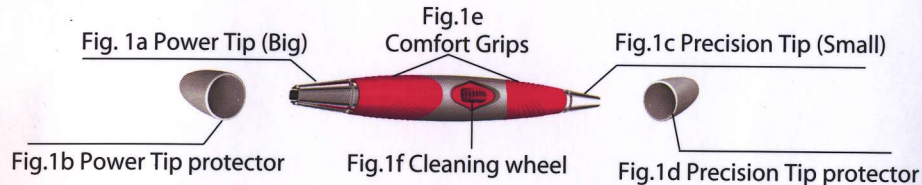


## Using your PediWow™ :

PediWow™ easily scrapes away the calluses, dry skin and corns and once they are removed the exfoliating tip will safely slide over the newly exposed soft and healthy skin.

**Before every use: The PediWow™ must be clean and free of all dead skin shavings.**

1. **Make sure that you fully understand all the instructions.**
2. Soak and wash in warm water the area to be treated thoroughly for about 5 -10 minutes (Fig. 2) making sure the dry and calloused skin is clean and soft. Lightly dry the area to be treated with a clean towel.

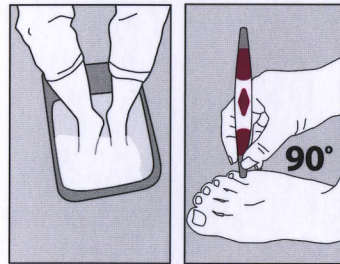

3. Remove the exfoliating plastic tip protectors (Fig. 1b) and (Fig. 1d). then choose the appropriate tip, either Power Tip for wide areas (Fig. 1a) or Precision Tip (Fig. 1c) for small or difficult to reach areas. Use smooth strokes on the area to be treated, securely holding the housing (Fig. 1) by the Comfort Grips (Fig. 1e). Holding the PediWow™ is as simple as holding a pen. It is very important to point the PediWow™ and the chosen tip at a straight angle (80 to 90 degrees) over the affected area (Fig. 3).


## ⚠ Important notice

\* **PediWow™** works perfectly on dry or wet skin


You may see different results according to the condition of the skin.

### On wet skin treatment

Dead skin shavings appear in the form of a wet paste or layers.

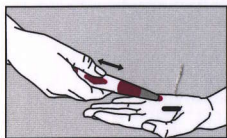


Wet paste

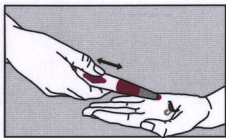

Layers

### On dry skin treatment

Dead skin shavings appear in the form of a thin dust.

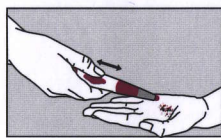

Thin dust

4. Now you may start applying light pressure and increasing it as needed. Work slowly, never rush the treatment process and stop frequently to check results. If the skin becomes irritated or if you feel pain, stop the treatment immediately. Once the healthy skin is revealed further treatment is no longer needed for that session. The callus shavings stored in the tip hole should be

emptied from time to time as needed by moving the cleaning wheel (Fig. 1f) in a backward motion with your thumb.


After treatment disinfect the skin then proceed to condition and moisturize it with a lotion to complete the **PediWow™** experience.

After using your **PediWow™** you should clean it using hot soapy water. If needed you may easily unscrew the stainless steel tips for deeper cleaning, finally rinse and dry it.


### Do and Don'ts



Straight smooth strokes=correct

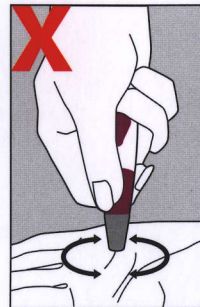

Twisting strokes=wrong

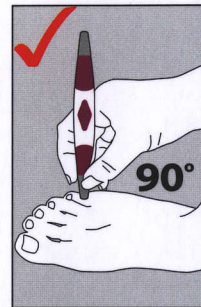

Vertical angle=correct



Slanted angle=wrong

PediWow™

Crystal nail file Fig. 2e ( 2 sides )

PediWow™

Nail buffer Fig. 2f ( 3sides ) 1 coarse side,  
1 fine side, 1 extrafine side

### PediWow kit is very easy to use

To use any of the following attachments: **Nail cleaner** Fig. 2a, **Cuticle pusher** Fig. 2b, **Ingrown toe nail file** Fig.2c and **Toe side nail file cleaner** Fig.2d, hold the **PediWow™** housing Fig.1, remove the Precision Cap Fig.1d and then unscrew the Precision Tip Fig.1c. Now you can easily screw in the chosen attachment. Apply pressure gently.

Better results are achieved if you soak the area in warm water for 5 to 10 minutes, then dry with a clean towel.

After using your **PediWow kit™** attachments you should clean them using hot soapy water. If needed you may unscrew them for deeper cleaning. Finally rinse, dry and store it.

Use the Crystal nail file Fig. 2e. to smooth out your nails edges. The Nail Buffer Fig. 2f. has 3 sides to buffer and polish your nails.

You may be able to carry the **PediWow™** and the

**PediWow™ kit** in the attractive Cosmetic bag Fig. 2g. and Tips case Fig. 2h.



PediWow kit carrying bag Fig. 2g



PediWow kit tips case Fig. 2h


## PediWow kit™ Instructions

### Using your PediWow kit™

Before using **PediWow kit™** make sure that you fully understand all the instructions.

The **PediWow kit™** attachments Fig. 2a, 2b, 2c and 2d are designed to be used just with **PediWow kit™** housing Fig.1

### PediWow™ housing Fig.1

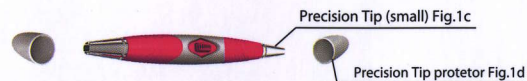

### PediWow kit™ components



Nail cleaner Fig. 2a



Cuticle pusher and cleaner Fig. 2b



Ingrown toe nail file Fig.2c



Toe side nail Fig. 2d

## Important Safety Instructions:

When using this device basic safety precautions should always be taken, including the following:

### Read all Instructions before using this device:

Important warnings:

- \* For intended use only. Intended for use in calloused only.
- \* Only use as described in this manual.
- \* Keep out of reach of children.
- \* Do not use on callus free smooth skin.
- \* Consult a physician before use if any skin condition exists.
- \* Consult a physician before use if you have any of the following but not limited to these symptoms: swelling, itchiness, pain or lack of sensitivity on the areas to be treated.
- \* Do not use on inflamed, irritated or infected skin.
- \* Do not use in open wounds.
- \* Do not use if you are diabetic, have poor blood circulation or have AIDS.
- \* Do not use if the metal becomes rusted, bent or damaged.
- \* Do not use if you are under chemotherapy treatment.
- \* Do not share this device with others.
- \* Discontinue use immediately if you experience any discomfort.